— ANTERIOR TRUNK SUPPORTS —

EN - USER INSTRUCTION MANUAL AND WARRANTY

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EN - INTRODUCTION

Thank you for purchasing a Ki Mobility product!

WARNING: Do not use this product without first reading this entire manual. Use without first reading this entire manual could result in serious injury to the user and/or attendant.

Ki Mobility recommends training by an Assistive Technology Professional (ATP) or other competent clinical or technical professional in the safe use of this product prior to use.

If you have any questions or concerns about any aspect of this product, this manual or the service provided, do not hesitate to contact your authorized supplier, an authorized representative or Ki Mobility. Ki Mobility can be reached at:

(1) 715-254-0991

Ki Mobility
5201 Woodward Drive
Stevens Point, WI 54481
U.S.A.

www.kimobility.com (Includes Dealer locator and access to printed information)

Or via our Authorized EU Representative:

Etac Supply Center AB
Långgatan 12
334 33 Anderstorp
Sweden

Or via our UK Responsible Person:

Etac Ltd
Unit 60, Hartlebury Trading
Estate, Hartlebury
Kidderminster, DY10 4JD

TABLE OF CONTENTS

EN - INTRODUCTION	1
A. Intended Use	1
B. Your Safety	1
C. Signal Words	1
D. Symbol Glossary	1
E. Check Your Skin Daily	2
F. General Warnings	2
G. Disposal	3
H. Re-Use	3
I. Storage	3
J. Maintenance	3
K. Cleaning	3
L. Sanitizing	3
M. Anterior Trunk Supports	4
KI MOBILITY LIMITED WARRANTY	5

AXIOM ANTERIOR TRUNK SUPPORTS

A. Intended Use

Ki Mobility Postural Support Devices are intended to provide postural support, positioning and improved sitting tolerance.

B. Your Safety

Any serious adverse events or injuries related to the use of this product, or its accessories must be reported immediately if the incident directly or indirectly leads to serious health decline or death of the user or other person. Report any serious events and/or injuries to Ki Mobility and, if required by local regulation, the competent health authority where the user and/or other person legally reside. Ki Mobility manufactures many different products that might meet your needs. The recognized best practice for selecting a product is to consult with an ATP or other competent technical professional and an experienced clinical professional such as physical therapist, occupational therapist, or physician. Final selection of the type of product, options and adjustments rests solely with you and your technical and clinical professionals. The options you choose, and the set-up and adjustment of the product have a direct impact on its performance, stability and its ability to meet your needs.

C. Signal Words

Within this manual you will find what are referred to as "Signal" words. These words are used to identify and convey the severity of varying hazards. Before using this product you, and each person who may assist you, should read this entire manual. Please note the Signal word and consider any notes, cautions or warnings. Make sure to follow all instructions and use your product safely. The Signal word refers to a hazard or unsafe practice that may cause severe injury or death to you or to other persons. The "Warnings" are in three main categories, as follows:

NOTE – Note indicates a potentially hazardous situation which, if not avoided, could result in a decline in the performance of the product or damage to the product.

CAUTION: Caution indicates a potentially hazardous situation which, if not avoided, could result in damage to the product and injury to you.

WARNING: Warning indicates a potentially hazardous situation which, if not avoided, could result in serious injury or death.

These signal words will be placed throughout the manual, where appropriate to highlight the hazardous situation.

D. Symbol Glossary

SYMBOL	DEFINITION
	Indicates the medical device manufacturer.
	Date of manufacture (YEAR/MM/DD).
	If presented on the medical device or packaging, it indicates the need for the user to consult the instructions for safety reasons such as cautions and warnings.
	Indicates the item is a medical device.
	Indicates the manufacturer's serial number so that a specific medical device can be identified.
	Indicates the authorized representative in the European Community.
	Indicates the manufacturer's declaration that the product meets the requirements of the applicable EC directives.
	Indicates the need for the user to consult the instructions for use.

D. Symbol Glossary (Continued)

SYMBOL	DEFINITION
	Indicates the need for the user to consult the listed website for instructions for use in an electronic format.
	Indicates a website where a user may obtain additional information about the medical product.
	Indicates a specified maximum weight limit (lbs/kg).
	Indicates a transit securement point.
	Indicates a potential pinch point.
	Indicates the entity importing the medical device into the locale.
	Indicates the manufacturer's batch code so that the batch or lot can be identified.
	Machine Wash, Hot
	Do Not Bleach
	Line dry or tumble, no heat
	Do Not Iron
	Do Not Dry Clean
	Indicates the manufacturer's declaration that the product meets the requirements of the applicable UKCA directives.

E. Check Your Skin Daily

These Axiom Seating and Positioning products are designed to improve pressure distribution and increase sitting tolerance. Whether the product provided will meet your individual needs cannot be assured in advance. The integrity of your skin can be affected by many aspects of your daily life and medical condition, including the use of this product. Be sure to follow any skin care regimens established by your clinician.

1. Users of this product should make sure their skin is inspected daily for changes. If you see changes, discontinue the use of this product immediately and consult with your clinician. A key clinical indicator of tissue breakdown is skin redness. If your skin develops redness, discontinue the use of this product immediately and consult your clinician. Be especially careful to inspect the areas related to the bony prominences of the pelvis such as your trochanters and ischial tuberosities.

F. General Warnings

WARNING: The postural support devices described in the user instruction manual were designed and are required to be fitted, installed, and adjusted for change exclusively by an Assistive Technology Professional (ATP or other competent clinical or technical professional). A clinical assessment should be done in advance of use to ensure safe and correct use. Improper selection or use could result in serious injury.

WARNING: Do not use this postural support device without first reading this entire manual. Use without first reading this entire manual could result in serious injury.

WARNING: Any caregivers must consult with the user's health care professional to learn proper adjustment, safety, and support methods best suited to the specific user.

WARNING: The installation of secondary supports on a wheelchair could affect the center of gravity of the wheelchair. Always review the instructions for use of your wheelchair to see and follow them to assure stability is not affected by the addition of any postural support devices. Changes in your center of gravity may affect your stability in your wheelchair, resulting in tipping over or falling from your wheelchair which may result in serious injury.

WARNING: Depending on your skin and other conditions related to your disability you may develop irritations or abrasions from prolonged contact with or rubbing against the surfaces of your wheelchair. Irritation and abrasions can lead to further complications from skin breakdown which can be serious.

Ki mobility recommends wearing clothing or other protected coverings over skin that is resting or rubbing against components or the frame of your wheelchair. If you have skin breakdown or other tissue integrity issues you should discontinue use and seek advice from a clinical professional trained in skin care such as a doctor, nurse, or a therapist. Always follow the advice of the clinical professionals caring for your skin and wear recommended dressings and/or coverings to protect your skin.

WARNING: Any body part that presses against a wheelchair frame or other wheelchair components, such as positioning devices, may create skin deformation that results in a high pressure point. High pressure points negatively affect the tissue under stress and can cause you to develop a pressure injury or skin breakdown. If your skin develops redness, or any other symptoms of a pressure injury or skin breakdown, discontinue the use of this product immediately and consult your clinician. Continuing to use the product after skin changes are present can result in serious injury or death.

WARNING: Failure to protect your Axiom secondary support(s) from extreme heat may cause your Axiom Product to ignite or emit harmful gases, resulting in serious injury or death.

WARNING: Your Axiom secondary support may contain small parts which can become disconnected or removed from the primary device and may pose a choking hazard, particularly to young children or individuals with cognitive impairments. Dispose of packaging and unused parts properly to prevent accidental ingestion or choking. Keep out of reach of children and vulnerable individuals and do not allow them handle or manipulate the device unsupervised. Failure to do so may result in accidental ingestion or choking which can result in serious injury or death.

G. Disposal

Product is made of materials which can be disposed of safely without special precautions.

Specific waste disposal or recycling regulations may be in force locally and these should be taken into consideration when disposal arrangements are made. This may include the cleaning or de-contamination of the component before disposal.

If recycling of materials is a requirement, please refer to the following list of general materials (which may or may not be present based on product type.):

- Aluminum – bracketry, clasps, clamps, hardware
- Steel – Fasteners, adjustment rods, extension rods
- Plastic – plugs, clips
- Upholstery - Woven polyester fabric and polyurethane foam.

If recycling of materials is desired or a requirement, such as per Directive 2012/19/EU, WEEE (Waste Electrical and Electronic Equipment) in Europe, both REAC and Mascot are compliant to WEEE and provide guidance on recycling and disposal. Proper dismantling is necessary to achieve a high proportion of reuse or recycle. Care to use proper tools and protective measures in handling are required to avoid contact with sharp edges. Plastics marked with recycling symbols and thermoplastic parts can be recycled. Aluminum, copper, and steel parts, including fasteners, may be recycled. Cables, batteries, and PCB are to be disposed of separately per local requirements. Return your product and accessories to your authorized supplier for proper disposal if you are not comfortable or able to dispose of properly.

H. Re-Use

In general, complete disinfection cannot be guaranteed on seams. We recommend that you dispose of upholstered components such as headrest pads and sewn support systems to avoid micro-bacterial contamination with active agents according to your local infection protection law.

I. Storage

Products should be stored in a cool dry place with relatively low humidity (<75%) out of direct sunlight and excessive dust for long periods of time. Products should not be stored outside.

J. Maintenance

WARNING: If you find any abnormalities, such as the foam breaking down or damage to the hardware or fasteners, discontinue the use of this product immediately and consult with your Authorized Ki Mobility Supplier. Failure to do so can affect your positioning and stability resulting in the potential for serious injury.

Regular cleaning and maintenance may help extend the life of your Postural Support Device. You should inspect your product monthly and during regular cleaning.

1. Check Cover for Tears and Excessive Wear

- Remove removeable covers and check for abnormalities (inside and out).
- Inspect the cover externally of pads with non-removeable covers.

iii. Check the foam, gel or other padding for breakdown or degradation.

1. Tears deeper than .25" (6.35mm)

2. Loose granules of foam

3. Permanent changes to the unloaded foam shape, such as flattening.

- Check the hardware for any looseness, signs of stripping, bending or cracking. Loose fasteners should be tightened to ensure product is securely in place. Consult your authorized dealer if the issue continues or if you have questions on appropriate retightening.
- If you find the fasteners or other hardware components are damaged and cannot be tightened, you should contact your Authorized Ki Mobility Supplier immediately to replace them and discontinue use of the product.
- If you find bent, cracked or otherwise damaged, such as corrosion, hardware components you should contact your Authorized Ki Mobility Supplier immediately to replace them and discontinue use of the product.

K. Cleaning

1. Removeable Product Covers

NOTE: Drying with heat may cause shrinking that affects the cover's fit and function.

- Remove the cover from product by detaching the Velcro® or unzipping the zipper on the cover.
- Remove the cover and turn the cover inside out.
- Zip the cover closed to avoid snagging.
- Machine wash in hot water (60°C).
- Line dry or tumble, no heat.

NOTE: Do not dry clean covers or use industrial washer or dryer to clean covers. Do not steam autoclave. Do not bleach. Laundry care symbol definitions are shown in image below.

Care & Use Instructions	
	Machine Wash, Hot
	Do Not Bleach
	Line dry or tumble, no heat
	Do Not Iron
	Do Not Dry Clean

2. Positioning Components

- Positioning components can be wiped clean with soap and water.
- DO NOT immerse positioning components underwater.

L. Sanitizing

Sanitizing step should be done in addition to regular cleaning:


- Remove and launder cover according to steps 1a., 1b. and 1c. above.
- Machine wash cover only in hot water (60°).
- Drip and tumble dry on No Heat.
- Disinfect all surfaces (cover, foam and back assembly) with over-the-counter disinfecting sanitizer of at least 70% alcohol or wipes. Do not soak or allow pooling of disinfectant solutions.

L. Sanitizing (Continued)

- e. Allow sanitizer to remain on surface for at least 15 minutes and remove with aseptic cloth if necessary.

M. Anterior Trunk Supports

Product Identification:

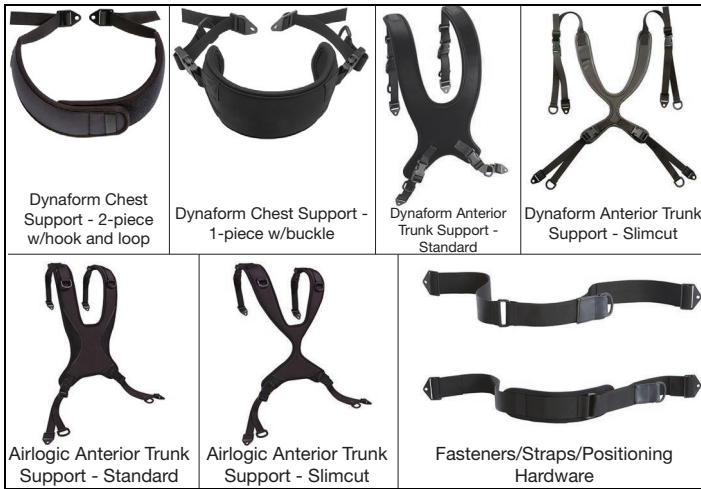⚠️ WARNING ⚠️

WARNING: Ensure restraint release mechanisms will not be activated by wheelchair components during a crash. Activation of restraint release mechanisms during a crash may result in serious injury or death.

- ⚠️ Use of headrests, lateral supports or other positioning accessories should not be used, or relied on as an occupant restraint. These belts should be positioned so that they don't interfere with the proper positioning of crash-worthy belt restraints and should not be relied on for occupant protection in crash situations unless the postural belt has been designed to comply with, and perform to, requirements 4.6, 5.2, 5.3 and 6.1 of RESNA WC-4.
- ⚠️ This secondary support is NOT intended as an occupant restraint in a moving vehicle. This component is for positioning use only and should not be relied upon as a transit safety restraint.
- ⚠️ Improper use of positioning belts and secondary supports can cause the user to slide underneath the positioning belt in the wheeled product. If this should happen, chest compression or suffocation due to pressure from the belt can occur. If this occurs, the user's breathing may be hampered causing serious injury or death.

Improper use of secondary positioning supports could cause body parts to become entangled which can restrict mobility and movement, which could result in serious injury or death. Secondary positioning supports can reduce the risk of falling from your wheeled product. Secondary positioning supports such as positioning belts, pelvic positioning belts and anterior trunk supports and straps, are designed to assist, retain and support proper positioning and posture in the wheeled product.

- Ensure the user does not slide underneath the secondary positioning support.
- The secondary positioning supports should have a snug fit; tight enough to hold their position, but not so

tight as to restrict breathing. An open hand should be able to fit between the belt and the user.

- Proper use of cushions can contribute to pelvic stability and reduce sliding.
- NEVER use secondary positioning supports and positioning belts:

- a. As a restraint, unless you can remove the belt easily in an emergency, then it may be considered a restraint. A restraint requires a doctor's order, consult with your health care advisor for other options to help you with your posture.
- b. On a user who is unconscious or agitated.
- c. As an occupant restraint or safety belt in a motor vehicle. A positioning belt or harness is not designed to replace a vehicle or wheelchair mount occupant restraint. During a sudden stop, with the force of the stop, the user would be thrown forward. Depending on the force, this may cause the positioning belt or harness to fail.

Use Instructions

⚠️ WARNING ⚠️


WARNING: Use of this postural secondary support without assuring stable pelvic positioning can result in the user sliding down or otherwise repositioning relative to the support and creating a choking hazard or other airway obstruction which could result in serious injury or death.

- ⚠️ Improper positioning of the support could become a choking hazard or other airway obstruction which could result in serious injury or death.
- ⚠️ Never use an anterior trunk support without appropriately installed pelvic positioning device(s) that assures pelvic stability. Failure to do so can result in slipping or repositioning which could result in serious injury or death.
- ⚠️ Never leave a user unattended during use of positioning support who is cognitively impaired or unable to free themselves from the device. Doing so could result in unintended repositioning, creating a choking hazard or other airway obstruction which could result in serious injury or death.
- ⚠️ Upper point of attachment must be above the shoulder for proper installation. This should be checked prior to use. Not following this instruction could result in unintended repositioning, creating a choking hazard or other airway obstruction which could result in serious injury or death.
- ⚠️ Ensure there is not pressure on the axilla (Armpit) area of client. This pressure area could cause discomfort and signifies your secondary support is not positioned correctly.

To be able to mount any anterior trunk support you must have a wheelchair with a solid backrest set approximately at shoulder height and installed pelvic positioning belt. Both wheelchair back and pelvic positioning belt must have been fit and installed by an authorized Ki Mobility supplier.

M. Anterior Trunk Supports (Continued)


1. Fitting and installation is done with the client in the product.
2. Ensure pelvic positioning belt is installed according to manufacturer's instructions with the oversight of the clinical professional.
3. Place straps of trunk support over the backrest and body of trunk support on chest of user.
4. Position strap between neck and shoulder joint.
5. Top Straps
 - a. Adjust each shoulder strap for a comfortable fit according to the user's shoulder width and sitting position.
 - b. Attach to the solid backrest using the fasteners / strap positioners.
6. Lower Straps
 - a. Identify optimal position for the lower straps on the backrest without interfering with other positioning devices (lateral supports, etc.).
 - b. Install end mounting between the top of the pelvis and the seat then attach using fasteners supplied (See image below).

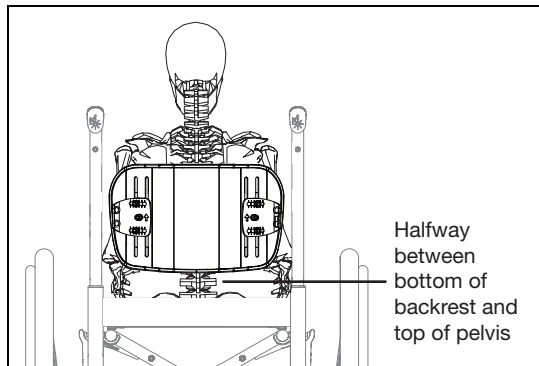7. Ensure system was safely installed by checking:
 - a. Buckles and straps are holding the support component in the desired clinically advised position.
 - b. User is sitting comfortably and in a way that reflects their normal seat posture without any part of the component obstructing airways, neck or head movement.

WARNING: Failure to read and heed these instructions could result in serious injury or death.

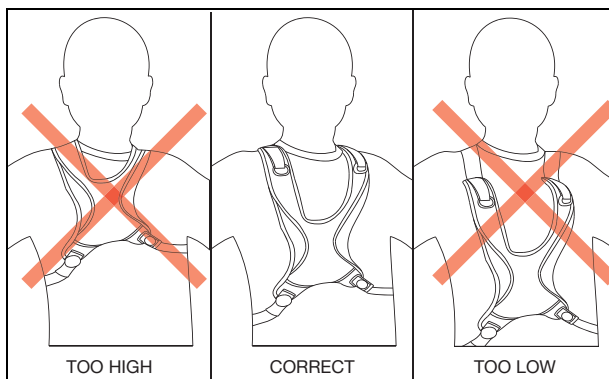KI MOBILITY LIMITED WARRANTY

The expected life of the product is 5 years.

A. For 24 Months

Ki Mobility warrants the Axiom Anterior Trunk Support to be free from defects in materials and workmanship for a period of 24 months. This limited warranty applies to:

1. Material and foam

For a period of 24 months (2 years) from the date of purchase, provided normal use, should a defect in materials or workmanship occur, Ki Mobility will, at its option, repair or replace it without charge. This warranty does not apply to punctures, tears, or burns.

B. For 6 Months

Ki Mobility warrants the Axiom Anterior Trunk Support to be free from defects in materials and workmanship for a period of 6 months. This limited warranty applies to:

1. Removable covers

For a period of 6 months from the date of purchase, provided normal use, should a defect in materials or workmanship occur, Ki Mobility will, at its option, repair or replace it without charge. This warranty does not apply to punctures, tears or burns.

C. Limitations

1. We do not warrant:
 - a. Damage from neglect, accident, misuse, or from improper installation or repair.
 - b. Products modified without Ki Mobility's express written consent.
 - c. Damage from exceeding the weight limit.
2. This warranty applies in North America only. Check with your supplier to find out if international warranties apply.
3. This warranty is not transferable and only applies to the first consumer purchase of this product through an authorized Ki Mobility supplier.

D. What We Will Do

Our sole liability is to repair or replace covered parts. This is the exclusive remedy for consequential damages.

E. What You Must Do

1. Obtain from us, while this warranty is in effect, prior approval for return or repair of covered parts.
2. Claims and repairs should be processed through the nearest authorized supplier. Except for express warranties made herein, all other warranties, including implied warranties of merchantability and warranties of fitness for a particular purpose are excluded.

F. Notice to Consumer

1. If allowed by law, this warranty is in place of any other warranty (written or oral, express or implied, including a warranty of merchantability, or fitness for a particular purpose).
2. This warranty gives you certain legal rights. You may also have other rights that vary from state to state.

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